

## MOZEC™ BALLOON TECHNICAL SPECIFICATIONS

Catheter System	:	Rapid Exchange (Rx)
Balloon Material	:	Nylon (Semi-Compliant)
Usable Catheter Length	:	142 cm
Max Guidewire Compatibility	:	0.014" (0.36 mm)
Guide Catheter Compatibility	:	5 F (min I. D. 0.056" / 1.42 mm)
Distal Shaft Coating	:	Biocompatible, Hydrophilic & Lubricious from distal Balloon base up to Rapid Exchange port
Balloon RO Markers	:	Platinum / Iridium 1 for Ø 1.25 to 1.50 mm, 2 for Ø 2.00 to 4.50 mm
Tip Length	:	5.00 mm for Ø 1.25 to 2.00 mm 3.50 mm for Ø 2.25 to 4.50 mm
Crossing Profile	:	0.60 to 1.50 mm
Proximal Shaft Diameter	:	1.98 F & 2.13F for length 45 & 49mm
Distal Shaft Diameter	:	2.40 F for Ø 1.25 to 2.00 mm 2.70 F for Ø 2.25 to 4.50 mm (Refer IFU for more details)
Markers on Proximal Shaft	:	Brachial and Femoral markers 90 cm and 100 cm from Distal Tip
Nominal Pressure (NP)	:	7 atm for all Diameters
Rated Burst Pressure (RBP)	:	16 atm Ø 1.25 to 4.00 mm 14 atm Ø 4.50 mm and for all Ø with length >41 mm (Refer IFU for more details)

## MOZEC™ ORDERING INFORMATION

Diameter\ Length	6 mm	9 mm	12 mm	14 mm	15 mm	17 mm	20 mm
1.25 mm	MOZ12506	MOZ12509	MOZ12512	-	MOZ12515	-	-
1.50 mm	-	MOZ15009	MOZ15012	-	MOZ15015	-	-
2.00 mm	-	MOZ20009	MOZ20012	MOZ20014	MOZ20015	MOZ20017	MOZ20020
2.25 mm	-	MOZ22509	MOZ22512	MOZ22514	-	MOZ22517	MOZ22520
2.50 mm	-	MOZ25009	MOZ25012	MOZ25014	-	MOZ25017	MOZ25020
2.75 mm	-	MOZ27509	MOZ27512	MOZ27514	-	MOZ27517	MOZ27520
3.00 mm	-	MOZ30009	MOZ30012	MOZ30014	-	MOZ30017	MOZ30020
3.25 mm	-	MOZ32509	MOZ32512	MOZ32514	-	MOZ32517	MOZ32520
3.50 mm	-	MOZ35009	MOZ35012	MOZ35014	-	MOZ35017	MOZ35020
4.00 mm	-	MOZ40009	MOZ40012	MOZ40014	-	MOZ40017	MOZ40020
4.50 mm	-	MOZ45009	MOZ45012	MOZ45014	-	MOZ45017	MOZ45020

Diameter\ Length	25 mm	30 mm	33 mm	38 mm	41 mm	45 mm	49 mm
2.00 mm	MOZ20025	MOZ20030	MOZ20033	MOZ20038	MOZ20041	MOZ20045	MOZ20049
2.25 mm	MOZ22525	MOZ22530	MOZ22533	MOZ22538	MOZ22541	MOZ22545	MOZ22549
2.50 mm	MOZ25025	MOZ25030	MOZ25033	MOZ25038	MOZ25041	MOZ25045	MOZ25049
2.75 mm	MOZ27525	MOZ27530	MOZ27533	MOZ27538	MOZ27541	MOZ27545	MOZ27549
3.00 mm	MOZ30025	MOZ30030	MOZ30033	MOZ30038	MOZ30041	MOZ30045	MOZ30049
3.25 mm	MOZ32525	MOZ32530	MOZ32533	MOZ32538	MOZ32541	MOZ32545	MOZ32549
3.50 mm	MOZ35025	MOZ35030	MOZ35033	MOZ35038	MOZ35041	MOZ35045	MOZ35049
4.00 mm	MOZ40025	MOZ40030	MOZ40033	MOZ40038	MOZ40041	MOZ40045	MOZ40049
4.50 mm	MOZ45025	MOZ45030	MOZ45033	MOZ45038	MOZ45041	MOZ45045	MOZ45049

Contact your country local Meril sales representative for availability of sizes highlighted in blue

Mozec™ Diameter 3.25mm and length 12mm is not CE approved  
Mozec™ NC is USFDA cleared (Refer IFU/Labels for approved sizes)

Mozec™ is registered trademark of Meril Life Sciences Pvt. Ltd.

Meril

More to Life

**Manufacturer:**  
**Meril Life Sciences Pvt. Ltd.**  
Survey No. 135/139, Bilakhia House,  
Muktanand Marg, Chala, Vapi - 396 191.  
Gujarat, India.  
T +91 260 305 2100  
F +91 260 305 2125

**Subsidiary companies:**  
**Meril Life Sciences Pvt. Ltd.**  
301, A-Wing, Business Square,  
Chakala, Andheri Kurla Road,  
Andheri East, Mumbai 400 093  
T +91 22 39350700  
F +91 22 39350777

**Meril, Inc.**  
2436 Emrick Boulevard,  
Bethlehem, PA 18020  
T +610 500 2080  
F +610 317 1672

**Meril South America**  
Doc Med LTDA  
Al. dos Tupiniquins,  
1079 - Cep: 04077-003 - Moema,  
Sao Paulo, Brazil.  
T +55 11 3624 5935  
F +55 11 3624 5936

**Meril GmbH.**  
Bornheimer Strasse 135-137,  
D-53119 Bonn,  
Germany.  
T +49 228 7100 4000  
F +49 228 7100 4001

**Meril Tibbi Cihazlar**  
Meril Tibbi Cihazlar İmalat ve Ticaret A.Ş.  
İçerenköy Mah. Çetinkaya Sok.  
Prestij Plaza No:28  
Kat:4 Ataşehir, 34752  
İstanbul / Turkey  
T +90 216 641 44 24  
F +90 216 641 44 25

**Meril China Co. Ltd.**  
2301b,23f, Lixin Plaza,  
no 90, South Hubin Road,  
Xiamen, China  
T 0086-592-5368505  
F 0086-592-5368519

**Meril SA Pty. Ltd.**  
102, 104, S101 and S102,  
Boulevard West Office Park,  
142 Western Service Road,  
Erf 813 Woodmead Extension 17  
Sandton, Johannesburg – 2191  
South Africa  
T +27 11 465-2049  
F +27 86 471 7941

**Meril Medical LLC.**  
Nauchnyi Proezd 19,  
Moscow, Russia – 117 246.  
Office - +7 495 772 7643

**EU representative.**  
Obelis S.A.  
Bd, General Wahis 53,  
1030, Brussels, Belgium.  
T +32 2 732 5954  
F +32 2 732 6003  
E mail@obelis.net

E askinfo@merillife.com  
W www.merillife.com

CE  
2460

MOZ/BROCHURE/009/20200201/GLOBAL

MOZEC™  
Rx PTCA Balloon Dilatation Catheter

Meril

MOZEC™

Rx PTCA Balloon Dilatation Catheter

FDA  
APPROVED


# MOZEC<sup>TM</sup>

Rx PTCA Balloon Dilatation Catheter

Handicraft for *Angiocraft*.


Mozec combines the competitive features of an ideal work-horse balloon dilatation catheter for a satisfying day-to-day cath-lab experience.

Hump free fusion of distal end for smooth transition

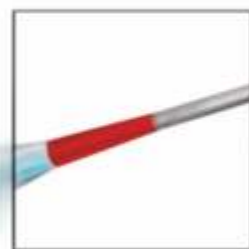


← Low tip entry profile & crossing profile

← Short and well matched balloon shoulders

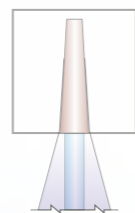

← Elongated 5 mm tip for Ø - 1.25, 1.50 & 2.00 mm for distal stability & guidewire support

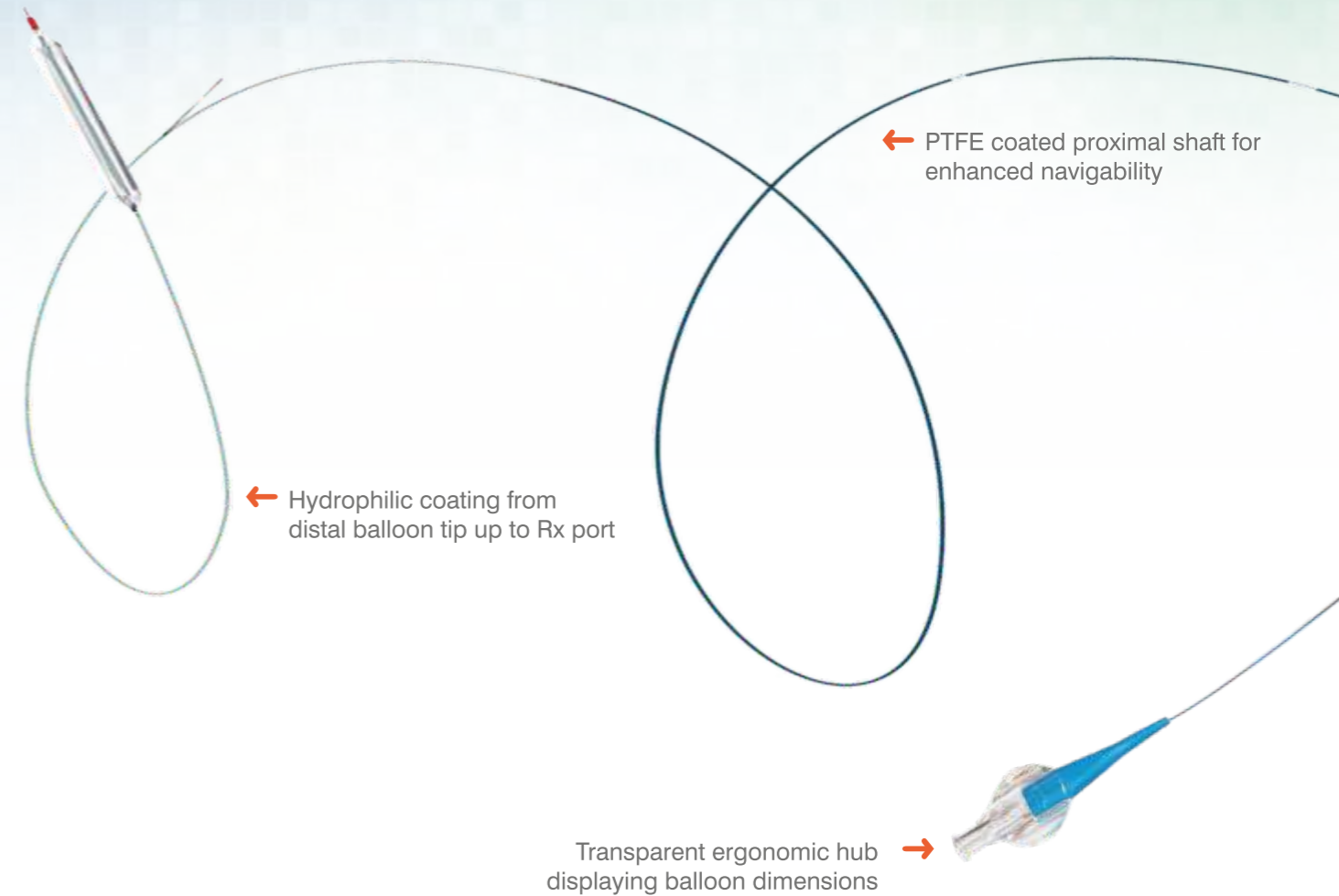

## Compliance Chart

Pressure (ATM)	1.25 mm	1.50 mm	2.00 mm	2.25 mm	2.50 mm	2.75 mm	3.00 mm	3.25 mm	3.50 mm	4.00 mm	4.50 mm
4	1.16	1.40	1.90	2.11	2.37	2.59	2.85	3.09	3.36	3.88	4.30
5	1.19	1.44	1.95	2.15	2.43	2.64	2.91	3.15	3.41	3.94	4.35
6	1.21	1.48	1.98	2.19	2.46	2.69	2.96	3.20	3.46	3.98	4.42
7	1.25	1.50	2.00	2.25	2.50	2.75	3.00	3.25	3.50	4.00	4.50
8	1.26	1.53	2.04	2.29	2.53	2.77	3.03	3.30	3.53	4.05	4.54
10	1.30	1.58	2.08	2.35	2.58	2.82	3.08	3.37	3.59	4.11	4.60
12	1.33	1.63	2.12	2.39	2.63	2.87	3.12	3.43	3.63	4.15	4.66
14	1.36	1.67	2.16	2.41	2.67	2.91	3.16	3.49	3.67	4.20	4.72
16	1.38	1.72	2.20	2.43	2.71	2.95	3.20	3.54	3.72	4.25	

● NP ● RBP